
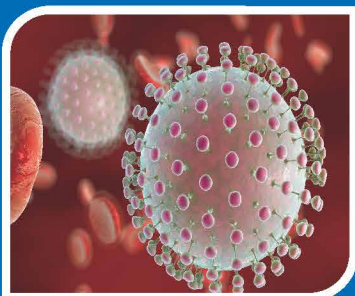


LumiQuick Diagnostics, Inc

In Vitro Diagnostics, Research,
Veterinary, Environmental, and more

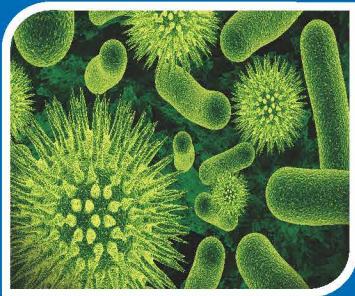

LumiQuick
DIAGNOSTICS, INC.

QuickProfile™
Quick Status®
Quicknostics®
EnviroSafe™

Drugs of Abuse
Infectious Diseases
Veterinary
Bio-threat Agents

Hormones
Cancer Markers
Cardiac Markers
ELISA

FDA cGMP ISO13485 CE



Company Introduction

LumiQuick Diagnostics Inc., located in the heart of Silicon Valley, California USA, develops, manufactures and markets high quality point of care tests and other immunoassay kits for the world-wide In vitro diagnostic market.

Through years of continuous and diligent efforts, we are proud to offer extensive product lines of rapid tests, ELISA, CLIA, Real Time PCR and reagent components. These product lines provide the immunoassays in various formats to detect infectious diseases, cardiac markers, cancer markers, drugs of abuse, fertility hormones and other disease markers.

Commitment with Quality

LumiQuick's products are manufactured per current Good Manufacturing Practices (cGMP) and ISO 13485.

Our mission is to provide the high quality products with affordable cost to help people around the world to fight human diseases and illicit substance abuse. We are committed to providing exceptional quality, customer service and cost saving to every customer.

We appreciate very much your interests in our products and look forward to be your partner for better healthcare.

QuickProfile™ Infectious Disease Rapid Test

Tropical and Hemorrhagic Fever Diseases



Cat. No.	Product Name	Qty/Kit
71019	Dengue IgG/IgM Combo Test Card	25
71066	Dengue NS1 Antigen Test Card	25
71087	Dengue IgG/IgM & NS1 Duo Test Card	25
71028	S. typhi Antigen Test Card	25
71054	S. typhi / paratyphi Antigen Test Card	25
71055	S. typhi IgG/IgM Duo Test Card	25
71055B	S. typhi IgG/IgM Combo Test Card	25
71049	Malaria pf Antigen Test Card	25
71050	Malaria pf/pv Antigen Test Card	25
71053	Malaria pan Antigen Test Card	25
71063	Malaria pf/pan Antigen Test Card	25
71017	Rickettsia IgG/IgM Combo Test Card	25
71018	Hantavirus IgG/IgM Combo Test Card	25
71031	Chikungunya IgG/IgM Combo Test Card	25
71078	Chagas Antibody Test Card	25
71090	Zika IgG/IgM Combo Test Card	25
71093	Zika NSI Test Card	25
71094	Zika IgG/IgM & NSI Duo Test	25

Respiratory Diseases

Cat. No.	Product Name	Qty/Kit
71034	Legionella Pneumonia Urinary Antigen Test Card	25
71048	Strep A Test Card	25
71062	Tuberculosis Antibody Test Card	25
71065	Influenza A (Flu A) Test Card	25
71067	Influenza A&B (Flu A&B) Combo Test Card	25

Gastro-Intestinal Diseases

Cat. No.	Product Name	Qty/Kit
71020	H. pylori Stool Antigen Test Card	25
71024	H. pylori Antibody Test Card	25
71029	Rotavirus Stool Antigen Test Card	25
71032	Adenovirus Stool Antigen Test Card	25
71033	Adeno/Rota Virus Stool Antigen Combo Test Card	25

QuickProfile™ Infectious Disease Rapid Test

Human Immunodeficiency Virus (HIV)


Cat. No.	Product Name	Qty/Kit
71001	HIV I&II Test Strip	50
71002	HIV I&II Test Card	25
71068	HIV I&II Tri-line Test Card	25



Hepatitis Virus

Cat. No.	Product Name	Qty/Kit
71003	HBsAg Test Strip	50
71004	HBsAg Test Card	25
71005	HBsAb Test Strip	50
71006	HBsAb Test Card	25
71007	HBcAb Test Strip	50
71008	HBcAb Test Card	25
71009	HBeAg Test Strip	50
71010	HBeAg Test Card	25
71011	HBeAb Test Strip	50
71012	HBeAb Test Card	25
71013-1	HBV-3 Panel Test Card (HBsAg/HBeAg/HBcAb)	25
71013-2	HBV-3 Panel Test Card (HBsAg/HBsAb/HBcAb)	25
71014	HBV-5 Panel Test Card (HBsAg/HBsAb/HBcAb/HBeAg/HBeAb)	25
71069-1	HBV-2 Panel Test Card (HBsAg/HBsAb)	25
71027	HCV Antibody Test Strip	50
71030	HCV Antibody Test card	25
71056	HAV IgG/IgM Test Card	25
71083	HAV IgM Test Card	25

*Customized combination of Panel Test Cards are available.

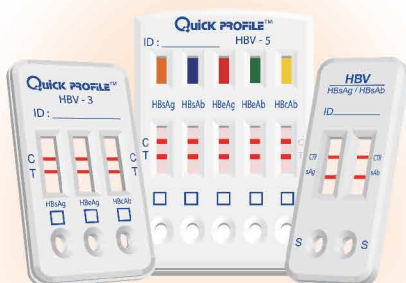

Sexually Transmitted Diseases (STD)

Cat. No.	Product Name	Qty/Kit
71015	Syphilis Antibody Test Strip	50
71016	Syphilis Antibody Test Card	25
71045	Chlamydia Test Card	25

QuickProfile™ Drugs of Abuse Rapid Test

Single Test Device

Cat. No.	Product Name	Qty/Kit
74001	Saliva Alcohol Test Strip	50
74049	Urine Alcohol Test Strip	50
74003 / 74002	Tramadol Test strip / Card	50 / 25
74013 / 74014	Amphetamine Test Strip / Card	50 / 25
74015 / 74016	Barbiturate Test Strip / Card	50 / 25
74017 / 74018	Benzodiazepine Test Strip / Card	50 / 25
74019 / 74020	Cocaine Test Strip / Card	50 / 25
74021 / 74022	EDDP Test Strip / Card	50 / 25
74023 / 74024	MDMA (Ecstasy) Test Strip / Card	50 / 25
74025 / 74026	Methadone Test Strip / Card	50 / 25
74027 / 74028	Methamphetamine Test Strip / Card	50 / 25
74029 / 74030	Opiate Test Strip / Card (300 ng/ml)	50 / 25
74031 / 74032	Opiate Test Strip / Card (2000 ng/ml)	50 / 25
74033 / 74034	Phencyclidine Test Strip / Card	50 / 25
74035 / 74036	Marijuana (THC) Test Strip / Card	50 / 25
74037 / 74038	Tricyclic Anti-Depressant (TCA) Test Strip / Card	50 / 25
74039 / 74040	Ketamine Test Strip / Card	50 / 25
74041 / 74042	Buprenorphine Test Strip / Card	50 / 25
74043 / 74044	Oxycodone Test Strip / Card	50 / 25
74051 / 74052	Propoxyphene Test Strip / Card	50 / 25
74055 / 74056	Methylphenidate (Ritalin) Test Strip / Card	50 / 25
74057 / 74058	Fentanyl Test Strip / Card	50 / 25
74059 / 74060	Clonazepam Test Strip / Card	50 / 25
74061 / 74062	Cotinine Test Strip / Card	50 / 25
74063 / 74064	Methaqualone Test Strip / Card	50 / 25
74065 / 74066	Synthetic Marijuana (K2) Test Strip / Card	50 / 25
74075 / 74076	Ethyl Glucuronide (ETG) Test Strip / Card	50 / 25
74077 / 74078	Synthetic Marijuana AB-Pinaca Test Strip/ Card	50 / 25
74079 / 74080	Synthetic Marijuana UR-144 Test Strip / Card	50 / 25
74081 / 74082	6-MAM (Heroin) Test Strip / Card	50 / 25



QuickProfile™ Drugs of Abuse Rapid Test



Multiple Drug Panel Test

Cat. No.	Product Name	Qty/Kit
74004 / 74004-TC	DOA-2 Panel, Dip Device / Cassette Device	25
74005 / 74005-TC	DOA-3 Panel, Dip Device / Cassette Device	25
74006 / 74006-TC	DOA-4 Panel, Dip Device / Cassette Device	25
74007 / 74007-TC	DOA-5 Panel, Dip Device / Cassette Device	25
74008 / 74008-TC	DOA-6 Panel, Dip Device / Cassette Device	25
74009 / 74009-TC	DOA-7 Panel, Dip Device / Cassette Device	25
74010 / 74010-TC	DOA-8 Panel, Dip Device / Cassette Device	25
74011 / 74011-TC	DOA-9 Panel, Dip Device / Cassette Device	25
74012 / 74012-TC	DOA-10 Panel, Dip Device / Cassette Device	25
74053 / 74053-TC	DOA-11 Panel, Dip Device / Cassette Device	25
74054 / 74054-TC	DOA-12 Panel Dip Device / Cassette Device	25
74045-A	Methamphetamine / THC Test Card	25
74045-B	Methamphetamine / Opiates Test Card	25
74046	DOA Panel Test Cup, split type, up to 14 drugs	25
74046-II	DOA Panel Test Cup, up to 14 drugs	25
74074	Saliva DOA Test Card, up to 9 drugs	25

QuickProfile™ Vitamin D

Rapid Test

Cat. No.	Product Name	Qty/Kit
83001	Vitamin D Test	25

ELISA

Cat. No.	Product Name	Qty/Kit
78065	Vitamin D EIA	96



QuickProfile™ Urine Strip

Cat. No.	Product Name	Qty/Kit
76001	Urine strip, up to 11 parameters	100
76002	7 Parameter Urine Adulteration Test Strip	25

Drugs of Abuse

Vitamin D

Urine Strips

QuickProfile™ Cardiac Markers Rapid Test

Cat. No.	Product Name	Qty/Kit
75003	Cardiac Panel Test Card, serum/plasma (TnI / CK-MB / Myoglobin)	25
75004	Cardiac Panel Test Card, serum/plasma/whole blood (TnI / CK-MB / Myoglobin)	25
75001	Troponin I Test Card, serum/plasma	25
75002	Troponin I Test Card, serum/plasma/whole blood	25
75005	Myoglobin Test Card, serum/plasma	25
75006	Myoglobin Test Card, serum/plasma/whole blood	25
75007	CK-MB Test Card, serum/plasma	25
75008	CK-MB Test Card, serum/plasma/whole blood	25

QuickProfile™ Cancer Markers Rapid Test

Cat. No.	Product Name	Qty/Kit
72001	Fecal Occult Blood Test Card	25
72002	Alpha Fetoprotein (AFP) Test Card	25
72003	Prostate Specific Antigen (PSA) Test Card	25
72009	CEA Test Card	25
72013	Transferrin Test Card	25
72014	Hemoglobin / Transferrin Test Card	25

QuickProfile™ Fertility Hormones Rapid Test

Cat. No.	Product Name	Qty/Kit
73001	HCG Midstream Test	25
73002	HCG Urine Test Card	25
73006	HCG Urine Test Strip	50
73005	HCG Serum/Urine Test Card	25
73007	HCG Serum/Urine Test Strip	50
73003	LH Midstream Test	25
73004	LH Urine Test Card	25
73008	LH Urine Test Strip	50



Cardiac Markers

Cancer Markers

Fertility Hormones

QuickProfile™ ELISA Kits

Cat. No.	Product Name	Qty/Kit
78001	HAV-IgM	96
78002	HAV-IgG	96
78003	HBsAg	96
78004	Anti-HBs (HBsAb)	96
78005	HBeAg	96
78006	Anti-HBe (HBeAb)	96
78007	Anti-HBc (HBcAb)	96
78008	Anti-HBc IgM (HBcIgM)	96
78009	HCV Antibody	96
78010	HIV Antibody	96
78011	HEV-IgM	96
78012	HEV-IgG	96
78013	HGV Antibody	96
78014	HTLV	96
78015	Treponema Pallidum (Syphilis)	96
78016	AFP	96
78017	CEA	96
78018	CA 125	96
78019	CA 15-3	96
78020	CA 19-9	96
78021	PSA	96
78022	Free PSA	96
78023	Free Beta-HCG	96
78024	Beta-2 Microglobulin	96
78025	TSH	96
78026	Ultra-Sensitive TSH	96
78027	T3	96
78028	Free T3	96
78029	T4	96
78030	Free T4	96
78031	Rubella IgG	96
78032	Rubella IgM	96
78033	Toxoplasma IgG	96
78034	Toxoplasma IgM	96
78035	CMV IgG	96
78036	CMV IgM	96
78037	HSV I IgG	96



Cat. No.	Product Name	Qty/Kit
78038	HSV I IgM	96
78039	HSV II IgG	96
78040	HSV II IgM	96
78041	H. pylori IgG, Qualitative	96
78042	H. pylori IgG, Quantitative	96
78043	H. pylori IgM	96
78044	HCG	96
78045	FSH	96
78046	LH	96
78047	Prolactin	96
78048	Visual HCG Pregnancy	96
78049	Ferritin	96
78050	HGH	96
78051	IgE	96
78052	Estradiol	96
78053	Progesterone	96
78054	Testosterone	96
78055	Troponin I	96
78056	Myoglobin	96
78057	High Sensitivity CRP	96
78058	Tuberculosis IgG	96
78059	HDV IgG	96
78060	HDV IgM	96
78061	HDV Antigen	96
78065	Vitamin D	96
78065A	Chlamydia IgA	96
78066G	Chlamydia IgG	96
78066M	Chlamydia IgM	96

QuickProfile™ Veterinary Rapid Test



Cat. No.	Product Name	Qty/Kit
9004	Canine Distemper Virus Antigen Test (CDV)	10
9005	Canine Parvovirus Antigen Test (CPV)	10
9006	Canine Coronavirus Antigen Test (CCV)	10
9007	Canine Rotavirus Antigen Test (CRV)	10
9008	Canine Adenovirus Type-II Antigen Test (CAV II)	10
9009	Canine Influenza Virus Antigen Test (CIV)	10
9010	Canine Parainfluenza Virus Antigen Test (CPIV)	10
9011	Canine Heartworm Antigen Test (CHW)	10
9012	Canine Leishmania Antibody Test (LSH Ab)	10
9013	Canine Ehrlichia canis Antibody Test	10
9014	Canine Brucella Antibody Test	10
9015	CPV/CCV Duo Test	10
9016	CCV/CRV Duo Test	10
9017	CDV/CAV II Duo Test	10
9018	CDV/ TOXO Duo Test	10
9019	CPV/CPIV Duo Test	10
9020	CPIV/CAV II Duo Test	10
9021	CDV/CPV Duo Test	10
9022	CPV/CDV/Giardia Panel Test	10
9023	CPV/CCV/CRV Panel Test	10
9024	CDV/CAV II/CPIV Panel Test	10
9025	CDV/CAV II/CIV Panel Test	10
9026	CDV/CAV I/CIV Panel Test	10
9027	Feline Panleukopenia Virus Antigen Test (FPV)	10
9028	Feline Infectious Peritonitis Antigen Test (FIPV)	10
9029	Feline Leukemia Virus Antigen Test (FeLV)	10
9030	Feline Immunodeficiency Virus Antibody Test (FIV)	10
9031	Feline Calicivirus Antigen Test (FCV)	10
9032	FIV/FeLV Duo Test	10
9033	Rabies Antigen Test	10
9034	Giardia Antigen Test	10
9035	Toxo Antigen Test	10
9036	Leptospira Antigen Test	10
9047	Porcine Circovirus Type II Antibody Test (PCV)	10
9048	Porcine Pseudorabies Virus Antibody Test (PPRV)	10
9049	Porcine Reproductive & Respiratory Syndrome Antibody Test	10
9050	Classical Swine Fever Virus Antibody Test	10

EnviroSafe™ Bio-Threat Agent Kits

Cat. No.	Product Name	Qty/Kit
9100	Anthrax Detection Kits	10
9101	Ricin Toxin Detection Kit	10
9102	Botulinum Toxin Detection Kit	10
9103	Y. pestis (Plague) Detection kits	10
9104	Staphylococcal Enterotoxin B (SEB) Detection Kit	10
9105	Tularemia Detection Kit	10
9106	Anthrax/Ricin/Botulinum/Y.pestis/SEB Panel Test	10
9107	Anthrax Training Kit	10
9108	Ricin Toxin Training Kit	10
9109	Botulinum Toxin Training Kit	10
9110	Staphylococcal Enterotoxin B (SEB) Training Kit	10
9111	Y.pestis (Plague) Training Kit	10
9113	Assorted Mix of Training Kit	10

Quicknostics® Real Time and Isothermal Amplification PCR

Cat. No.	Product Name	Qty/Kit
82001	HPV, high risk, isothermal real time	96/48
82002	HPV, high risk, isothermal real time with simultaneous genotyping HPV 16 and 18	96/48
82003	HPV, high risk, isothermal real time, genotyping all 15 high risk HPV	96/48
82004	HPV, high risk, real time PCR	96/48
82005	HPV, high risk, real time PCR with simultaneous genotyping HPV 16 and 18	96/48
82006	HPV, high risk, real time PCR, genotyping all 15 high risk HPV	96/48
82007	Chlamydia trachomatis, real time PCR	96/48
82008	Clostridium difficile, real time PCR	96/48
82009	Coxiella burnetii, real time PCR	96/48
82010	Mycoplasma pneumonia, real time PCR	96/48
82011	Neisseria gonorrhoeae, real time PCR	96/48
82012	Yersinia pestis, real time PCR	96/48
82013	Bacillus Anthracis, real time PCR	96/48
82014	Francisella tularensis, real time PCR	96/48
82015	Haemophilus, real time PCR	96/48
82016	Legionella, real time PCR	96/48
82017	Listeria, real time PCR	96/48
82018	Salmonella, real time PCR	96/48
82019	Staphylococcus, real time PCR	96/48
82020	Streptococcus, real time PCR	96/48
82021	Trichomonas vaginalis, real time PCR	96/48

QuickProfile OxiStress™ Tests



Cat. No.	Product Name	Qty/Kit
9201	Lactoferrin EIA Test Kit	10
9202	Aconitase Activity Colorimetric Assay Kit	10
9203	Urinary 8-epi-Prostaglandin F2a Assay Kit	10
9204	Dityrosine EIA Test Kit	10
9211	Reduced glutathione / oxidized glutathione Glutathione Colorimetric Assay Kit	10
9211MP	Reduced glutathione / oxidized glutathione Glutathione Colorimetric Assay Kit	10
9212	lutathione Reduced Colorimetric Assay Kit	10
9212MP	Glutathione reduced microplate assay kit	10
9213	Cellular Glutathione Peroxidase Colorimetric Assay Kit	10
9214	Glutathione Reductase Colorimetric Assay Kit	10
9215	High sensitivity Glutathione Reduced Colorimetric Assay Kit	10
9215MP	High sensitivity Glutathione Microplate Assay Kit	10
9216	Catalase Colorimetric Assay Kit	10
7217	Hydrogen Peroxide Colorimetric Assay Kit	10
9218	Cellular Glutathione Peroxidase Control	10
9219	Glutathione Reductase Control	10
9221	Hexanoyl-Lysine EIA kit	10
9222	Lipid Peroxidation Colorimetric Assay Kit	10
9222MP	Lipid Peroxidase Microplate assay kit	10
9223	Malondaldehyde Colorimetric Assay Kit	10
9223MP	High sensitivity malondaldehyde microplate assay kit	10
9224	2-Thiobarbituric Acid Reactive Substances Colorimetric Assay Kit	10
9224MP	2-Thiobarbituric Acid Reactive Substances Microplate Assay Kit	10
9225	Malondaldehyde microplate assay kit	10
9231	DNA Damage Quantitative Microplate Assay Kit	10
9232	Regular 8-OHdG-EIA Kit	10
9241	Antioxidant Potential Colorimetric Assay Kit	10
9261	Ultra sensitive Assay for Nitric Oxide Synthase	10
9262	Nitric Oxide Non-Enzymatic Assay Kit	10
9263	Colorimetric Nitric Oxide Assay Kit	10

Quickprofile™ Rapid Test

Adeno-Rotavirus



Quick Profile™ Adeno/Rota Combo Test is an in vitro qualitative immunochromatographic assay for the rapid detection of adenovirus and rotavirus antigens in human stool specimen. The test results are intended to aid in the diagnosis of adenovirus and rotavirus infection and to monitor the effectiveness of therapeutic treatment.

Specimen: Stool

Assay Reading Time: 10 minutes

Accuracy:

Adenovirus: Pos.=96.57%, Neg.=99.29%

Rotavirus: Pos.=97.27%, Neg.=99.28%

Storage: 4 to 30°C

Shelf Life: Shelf Life: 18 months from production

Chikungunya IgG/IgM



QuickProfile™ ChikV IgG/IgM Combo Test is a rapid immunochromatographic assay for the simultaneous detection of IgG and IgM antibodies to Chikungunya virus in human whole blood, serum or plasma. The assay is used as a screening test for Chikungunya viral infection.

Specimen: Whole blood/Serum/Plasma

Assay Reading Time: 15 minutes

Accuracy:

IgG/IgM: POS.=96.1%, Neg.=91.67%

Storage: 4 to 30°C

Shelf Life: 18 months from production

Dengue NS1 Ag & IgG/IgM Ab



QuickProfile™ Dengue Duo Panel is an in vitro qualitative immunochromatographic assay for the rapid detection of Dengue NS1 antigens and anti-dengue IgG/IgM antibodies in human blood, serum and plasma specimens simultaneously. The test results are intended to aid in the diagnosis of dengue infection.

Specimen: Whole blood/Serum/Plasma

Assay Reading Time: 20 minutes

Accuracy:

Dengue IgG: Pos.=91.43%, Neg.=98.67%

Dengue IgM: Pos.=95.40%, Neg.=98.85%

Dengue NS1: Pos.=93.6%, Neg.=100%

Storage: 4 to 30°C

Shelf Life: 18 months from production

Influenza A & B



QuickProfile™ Influenza A+B Test is an in vitro immunochromatographic assay for the qualitative detection of influenza type A and type B nucleoprotein antigens in nasopharyngeal swab, nasal swab and nasal aspirate samples. The identification is based on the monoclonal antibodies specific for the nucleoprotein of Influenza virus A and B.

Specimen: Whole blood/Serum/Plasma

Assay Reading Time: 10 minutes

Accuracy:

Flu A: Pos.=91.7%, Neg.=98.07%

Flu B: Pos.=91.7%, Neg.=98.07%

Storage: 4 to 30°C

Shelf Life: 18 months from production

Quickprofile™ Rapid Test

HBsAg



QuickProfile™ HBsAg Test is an immunochromatography assay for the qualitative detection of Hepatitis B virus surface antigen (HBsAg) in human serum, plasma or whole blood specimen.

Specimen: Whole blood/Serum/Plasma

Assay Reading Time: 20 minutes

Accuracy:

Pos.=99.6%, Neg.=99.6%

Storage: 4 to 30°C

Shelf Life: 18 months from production

H. pylori Ag



QuickProfile™ H. pylori Antigen Test is an in vitro qualitative immunochromatographic assay for the rapid detection of Helicobacter pylori antigens in human stool specimen. The test results are intended to aid in the diagnosis of H. pylori infection, to monitor the effectiveness of therapeutic treatment and to confirm the eradication of H. pylori in peptic ulcer patients.

Specimen: Stool

Assay Reading Time: 10 minutes

Accuracy: Pos.=99.03%, Neg.=97.95%

Storage: 4 to 30°C

Shelf Life: 18 months from production

Zika IgG/IgM



QuickProfile™ Zika Virus IgG/IgM Combo Test is a rapid immunochromatographic assay for the simultaneous detection of IgG and IgM antibodies to Zika virus in human whole blood, serum or plasma. The assay is intended for diagnostic use only.

Specimen: Whole blood/Serum/Plasma

Assay Reading Time: 15 minutes

Accuracy:

IgG/IgM: Pos.=92.31%, Neg.=90.91%

Storage: 4 to 30°C

Shelf Life: 18 months from production

Cardiac Markers Panel



Quick Profile™ Cardiac Panel test is an immunochromatography based one step in vitro test. It is designed for qualitative determination of cardiac troponin I (cTnI), CK-MB and Myoglobin in human serum, plasma or whole blood specimens as an aid in the diagnosis of myocardial infarction.

Specimen: Whole blood/Serum/Plasma

Assay Reading Time: 15 minutes

Accuracy:

Troponin I: Pos.=100% Neg.=100%

CK-MB: Pos.=100% Neg.=99.3%

Myoglobin: Pos.=100% Neg.=99.3%

Storage: 4 to 30°C

Shelf Life: 18 months from production

Quickprofile™ Rapid Test

HCG Pregnancy



Quick Profile™ HCG Pregnancy Urine Test is a chromatographic immunoassay for the early detection of human chorionic gonadotropin (HCG) in urine specimens. The test is intended for healthcare professional use.

Specimen: Urine

Assay Reading Time: 5 minutes

Sensitivity: 25mIU/mL

Accuracy: Pos.=100%, Neg.=99.5%

Storage: 4 to 30°C

Shelf Life: 24 months from production

Drugs of Abuse



QuickProfile™ DOA/Alcohol Panel Test is an immunochromatography based one step in vitro test. It is designed for qualitative determination of drug substances in human urine specimens. This assay may be used in the point of care setting.

Specimen: Urine

Assay Reading Time: 5 minutes

Storage: 4 to 30°C

Shelf Life: 18 months from production

Accuracy Agreement			
Component	Cut-off	POSITIVE	NEGATIVE
Amphetamine	1000 ng/mL	98.5%	100%
Barbiturate	300 ng/mL	100%	100%
Benzodiazepine	300 ng/mL	98%	100%
Buprenorphine	10 ng/mL	96%	100%
Cocaine	300 ng/mL	99%	99%
Clonazepam	300 ng/mL	100%	100%
Cotinine	100 ng/mL	100%	100%
EDDP	100 ng/mL	98%	100%
Fentanyl	10 ng/mL	100%	100%
K2	50 ng/mL	100%	100%
Ketamine	100 ng/mL	99%	100%
MDMA	500 ng/mL	96%	95%
Methadone	300 ng/mL	100%	100%
Methamphetamine	1000 ng/mL	98%	100%
Methylphenidate	300 ng/mL	100%	100%
Opiate	300 ng/mL	99%	99%
Opiate II	2000 ng/mL	94%	100%
Oxycodone	100 ng/mL	100%	95%
Phencyclidine	25 ng/mL	98%	95%
Propoxyphene	300 ng/mL	100%	100%
THC	50 ng/mL	100%	99%
Tramadol	200 ng/mL	95%	98%
TCA	1000 ng/mL	98%	95%



QuickProfile OxiStress™ Tests

LumiQuick Quick Profile OxiStress™ is designed to provide a quick and convenient tool for the measurement of oxidative stress in humans and animals, through the measurement of markers such as protein, enzyme, lipid & DNA.



To provide reference information on human health status.

Improve livestock's productivity and product quality.

Study of oxidative mechanism in bio-specimens to provide the clue for curing the diseases caused by oxidative stress.

Diseases Often Associated to Oxidative Stress:

ADHD, Cancer, Parkinson's disease, Lafora disease, Alzheimer's disease, Atherosclerosis, Heart Failure, Myocardial Infarction, Fragile X Syndrome, Sickle-cell Disease, Lichen planus, Vitiligo, Autism, Infection, Chronic fatigue syndrome and Depression.



LumiQuick
DIAGNOSTICS, INC.

**POCT
CLIA Waived
Forensic Use
Manufactured in the USA**

Quick Profile™ Drugs of Abuse Tests



**Test Cup - Test Card - Dip Test - Multi-panel
Drug Tests available in individual test or combo panels**

AMPHETAMINE
BARBITURATE
BENZODIAZEPINE
BUPRENORPHINE
COCAINE
METHADONE
METHAMPHETAMINE
OPIATES
OXYCODONE
PCP
THC
ADULTERATION STRIPS

FDA approved

FENTANYL
TRAMADOL
CLONAZEPAM/7-ACL
COTININE
EDDP
ETG
K2 (JWH, AB-PINACA, UR-144)
KETAMINE
6-MAM (HEROIN)
MDMA/ECTASY
METHYLPHENIDATE (RITALIN)
PROPOXYPHENE
TCA
ALCOHOL (URINE OR SALIVA)

Forensic Use Only



Test Operation

<p>1</p>	<p>1</p>	<p>1</p>	<p>1</p>
<p>2</p>	<p>2</p>	<p>2</p>	<p>2</p>
<p>3</p>	<p>3</p>	<p>3</p>	<p>3</p>
<p>4</p>	<p>4</p>	<p>4</p>	<p>4</p>



Test Operation

<p>1</p>
<p>2</p>
<p>3</p>
<p>4</p>

Test Result Appearance and Interpretation

<p>NEGATIVE</p>	<p>POSITIVE</p>	<p>INVALID</p>
------------------------	------------------------	-----------------------

Ordering Information


- Quick Profile™ Single Drug Strip or Card Test
- Quick Profile™ Multiple Drug Test Dip Device or Card Device
- Quick Profile™ Multiple Drug Test Cup
- Quick Check™ Adulteration Strip Test

www.quickprofileDOA.com

Follow us on social media:



2946 Scott Blvd, Santa Clara CA 95054 USA
 Tel.: +1(888)352-8394 - info@LumiQuick.com
www.LumiQuick.com

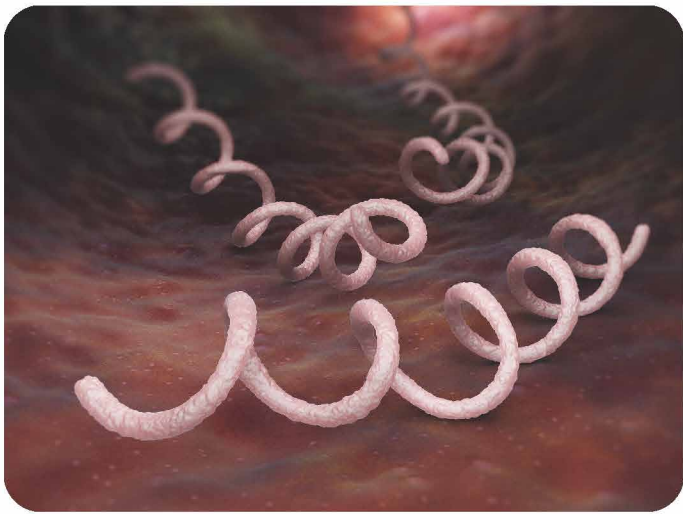

Quick Profile™ Syphilis Antibody Test

(Chromatographic Lateral Flow)

Quick Profile™ Syphilis Antibody Test is a rapid immuno-chromatographic assay for the qualitative determination of *Treponema pallidum* antibodies.

Specimen:
Human Whole Blood
Serum or Plasma.



Quick Profile™ Syphilis Antibody
One Step Test available as Test Card or Test Strip

Quick Profile™ Syphilis Antibody Test


Procedure Test Card

Bring the kit components to room temperature.

Transfer 150µL of sample using the pipette provided or a micropipette.

Hold the pipette in a vertical position over the sample well and deliver 3 drops (120-150µL).

Read the results at 15 min.

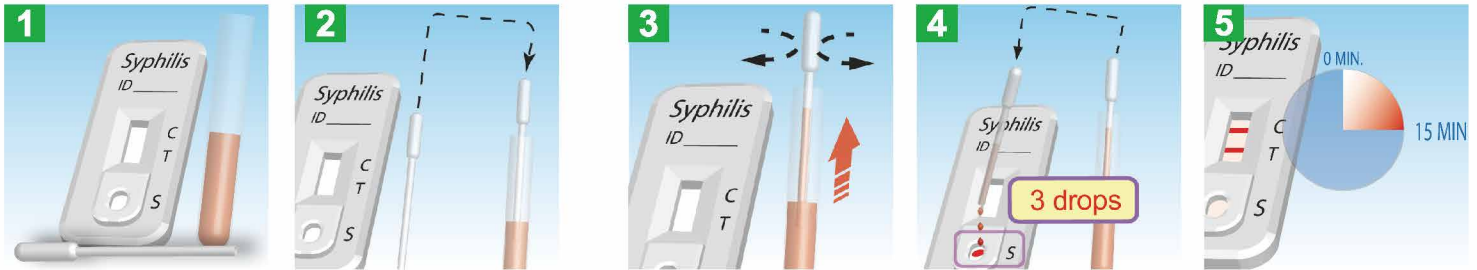


Procedure Test Strip

Bring the kit components to room temperature.

Dip the strip into the sample with the arrow marked end down.

Keep the strip in touch with the sample until it reaches the membrane. Remove the strip and place it on a flat surface.

Read the results at 15 min.

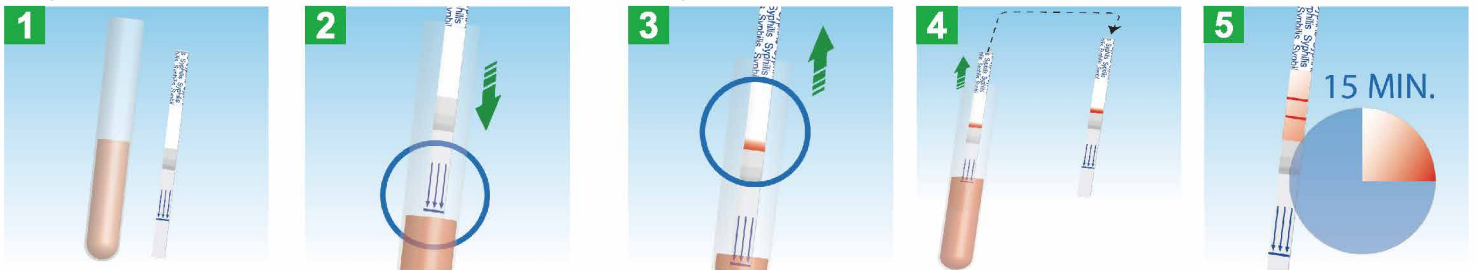

Test Result Appearance and Interpretation



Product Performance

REFERENCE	QuickProfile™ Syphilis Antibody Test	
	POSITIVE	NEGATIVE
TEST RESULTS		
POSITIVE SAMPLE	349	7
NEGATIVE SAMPLE	5	709
TOTAL SAMPLES	354	716
ACCURACY	98.6%	99%

Ordering Information

Catalog No	Description	Qty of test per Kit
71015	Quick Profile™ Syphilis Antibody Test Strip	50 tests
71016	Quick Profile™ Syphilis Antibody Test Card	25 tests



2946 Scott Blvd, Santa Clara CA 95054 USA
 Tel.: +1(408)855-0061 - info@LumiQuick.com
 www.LumiQuick.com

LumiQuick





DIAGNOSTICS, INC.

2946 Scott Blvd.
Santa Clara, California
95054 USA
www.LumiQuick.com



www.youtube.com/lumiquickinc

YouTube

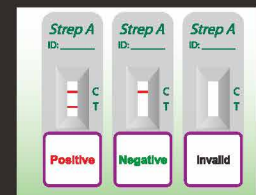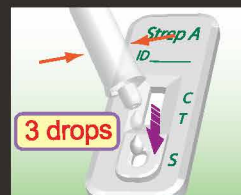

Questions? Please contact us at
info@lumiquick.com

Or, your local agent 

For Authorized Distributors Use